



DATAPAQ™ MULTIRACK

2D Barcode Scanner

For All Rack and Box Formats

FEATURES

- Capable of scanning Cryobox formats
- Capable of scanning SBS Rack formats
- Compatible with all 2D Datamatrix tubes
- 5 mask types allow a variety of rack and box combinations to be scanned
- DataPac™ software with free lifetime upgrades
- Cryoprotection option available

Easy to Set-Up and Easy to Use

The DataPac™ MultiRack 2D scanner offers easy, out-of-the-box set up and is delivered pre-calibrated and ready to scan.

This instrument can handle up to 4 SBS racks or 2 Cryoboxes in a single scan. It is supplied with all the features needed for robotic integration.

The MultiRack is a good option for those wanting flexibility; it can be used for 4 x SBS racks one day, and a single Cryobox the next.

SPECIFICATIONS

- Dimensions = 270mm x 475mm x 75mm (WxDxH)
- Total Weight = 3.1kg
- Power input = 100-240VAC
- Compatible with 24, 48, 96 position SBS racks
- Compatible with 81 and 100 position Cryoboxes
- Compatible with all commercially available datamatrix 2D coded tubes
- Export data to Excel, text, JSON, XML and image files
- Operating system on host PC = Windows 7, 8, 10
- Data output connector = USB
- Warranty period = 2 years
- Free DataPac™ software upgrades for life
- Software includes remote control capability for integration applications

INCLUDED WITH THE MULTIRACK

- DataPac™ MultiRack instrument
- Free choice of one mask
- Power cable
- DataPac™ Software
- USB lead
- User guide




4 x SBS mask is supplied with mirrors that allow the MultiRack to read the 1D rack codes



OPTION

CRYOPROTECTION

Prevents the scanner window from misting up when a cold rack is scanned. Our passive technology does not use heat or blown air

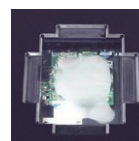

Ordering information

Code	Description	Quantity
ZTS-A4	DataPac MultiRack Scanner	1
ZTS-A4-80	DataPac MultiRack with Cryoprotection	1
ZTS-MSBCR	Mask for 1 x SBS Rack and 1 x Cryobox	1
ZTS-MSB4	Mask for 4 x SBS Rack	1
ZTS-MSB3	Mask for 3 x SBS Rack	1
ZTS-MCRYO2	Mask for 2 x Cryoboxes	1
ZTS-MCRYO1	Mask for 1 x Cryobox	1